


RespiTrainer®  
Advance

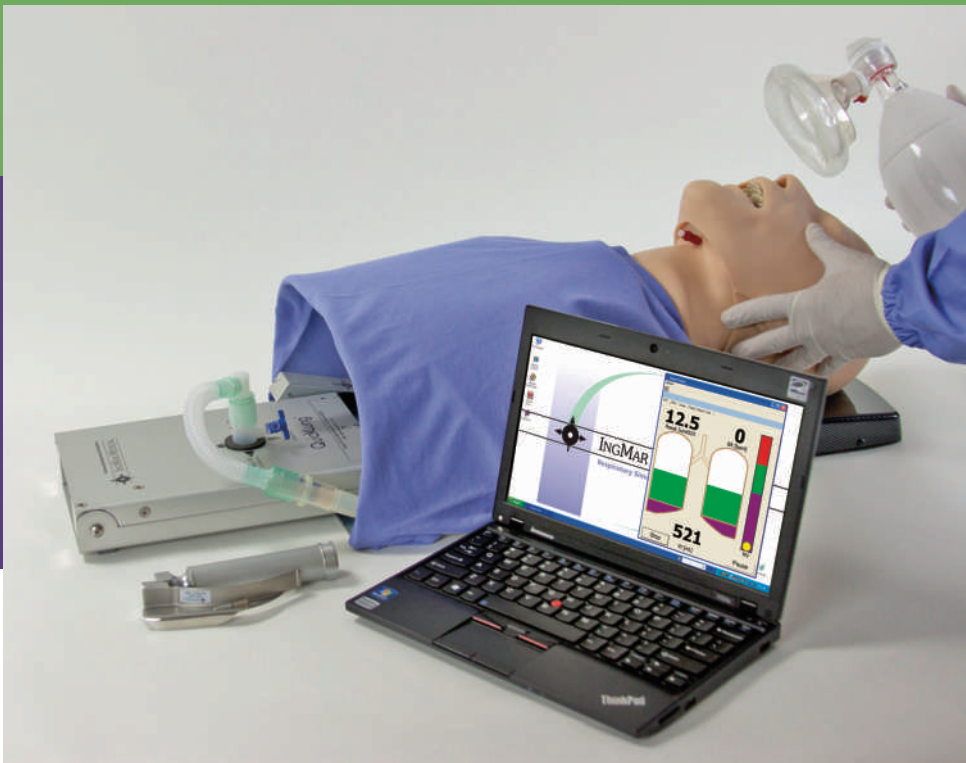

# Multi-skill Ventilation and Airway Management

If your students or staff were less than adequately trained in ventilation and endotracheal intubation, how would you know? RespiTrainer Advance puts you inside your patient's airway and lungs with instant feedback on delivered ventilation.

With RespiTrainer Advance, you get the essential features of a high-end respiratory manikin plus true ventilation data – at a fraction of the cost. Easy set-up for both instructor and student, so you can start training and testing without delay.


The adjustable lung enables hands-on experience with different patient lung conditions to prepare for a wide range of real life scenarios. Chest rise feature provides realistic visual feedback.

The RespiTrainer Advance training head offers unparalleled airway realism for LMA and ET-tube insertion training. Features include anatomically correct airway with visually accurate features, an inflatable tongue to increase difficulty of airway, break-out teeth, and "real-feel" skin covering for practicing proper mask placement.



**ETView** – See exactly what the student sees during intubation training. A tiny digital camera at the tip of the ET-tube transmits video images of RespiTrainer Advance's airways.

# RespiTrainer® – Feedback for Proper Training

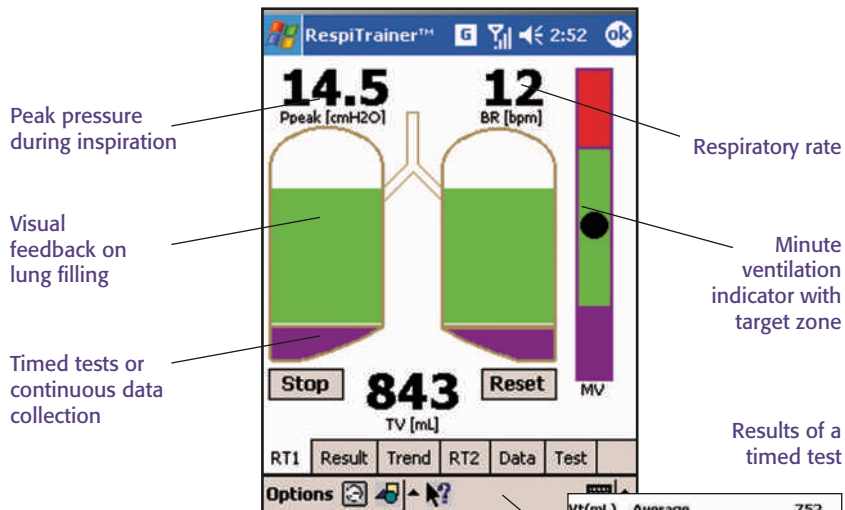

Breath by breath spreadsheet analysis of ventilation trial can be used for competency testing.

Student: Jane Doe  
One minute manual bag mask ventilation  
RespiTrainer Data File

Time Sec.	BREATH	BrRate	Vt	Peak Pressure
0.0	START 1			
4.2	1	14.3	290	6.0
8.1	2	15.4	313	6.4
12.6	3	13.3	423	8.6
16.9	4	14.0	561	11.5
21.4	5	13.3	524	10.7
26.6	6	11.5	544	11.1
31.9	7	11.3	580	11.8
36.3	8	13.6	479	9.8
41.2	9	12.2	562	11.4
46.0	10	12.5	523	10.7
51.0	11	12.0	497	10.2
55.6	12	13.0	511	10.4
60.0	END 1			

#### Results

	Br	Vt	Pp	MV
60.0 AVE	13.0	484	9.9	5665
<b>Thresholds</b>				
60.0 >MAX 1	0	0	0	
60.0 <MIN	0	2	0	0

## Related Products

In addition to **ETView®** (pictured on the front), the following options are also available:


### RespiTrainer Infant

Simulate your tiniest patients.

### Breather Option

Simulate the spontaneously breathing patient in a variety of modes and patterns – a capability previously available only on high-end breathing simulators.

© 2013 IngMar Medical. RespiTrainer is a trademark of IngMar Medical, Ltd. ETView is a trademark of ETView Ltd.

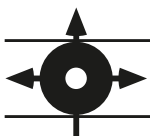

# INGMAR MEDICAL

Respiratory Simulation Specialists

IngMar Medical, Ltd. is ISO 9001:2008 certified.

5940 Baum Boulevard  
Pittsburgh, PA 15206 USA  
412.441.8228 Toll free 800.583.9910  
www.ingmarmed.com

© 2015 IngMar Medical, Ltd.